










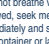


SBio ABD Combipack

Anti-A, Anti-B, Anti-D (Rho) (IgM+IgG) Reagent Combipack For Routine Blood Grouping And Typing

	Anti-A	Anti-B	Anti-D (Rho) IgM+IgG
REF	90220310		
Pack	3 x 10 ml		

 + 2°C Store at 2-8°C	 Manufacturer	 In vitro Diagnostic Medical Device	 Batch Number / Lot Number	 Expiry date	 Xn H302, H332 S23-S24/25
 Consult Instructions for use	 Date of Manufacture	 Catalogue Number	 This side up	 Description of the reagent	 Harmful if swallowed. Do not breathe vapour. If swallowed, seek medical advice immediately and show this container or label. Avoid release to the environment. Refer to special instructions.

INTENDED USE

SBio ABD Combipack consists of monoclonal blood grouping Anti-A, Anti-B, Anti-D (Rho) (IgM + IgG) reagents. These reagents can be used for *in vitro* determination of the presence or absence of A, B, D (Rho) antigen and D^u antigen, a variant of D (Rho) antigen, found on human red blood cells.

SUMMARY

Monoclonal antibodies are derived from hybridoma cell lines, created by fusing mouse antibody producing B lymphocytes with mouse myeloma cells or through EBV transformation. Each hybridoma cell line produces homogenous antibodies of only one immunoglobulin class, which are identical in their chemical structure and immunological activity. Since monoclonal antibodies only recognize one epitope, they generally have low cross reactivity with non-specific antigens. Their epitope specificity, limited cross reactivity, and long term yield make monoclonal antibodies attractive for use in many biological assays and application.

In blood transfusion the two most important blood group systems are the ABO system and the Rhesus system. The ABO system divides the blood groups into four major groups A, B, AB and O depending on the presence or absence of the A or B antigens on the red blood cells. Approximately 41% of the Caucasian population have the A Antigen, 9% have the B Antigen, 4% have both A and B antigens, while the remaining have neither A nor the B antigen.

The second major blood group system is the Rhesus (Rho) system. D is the main antigen of the Rhesus system. Human red blood cells are classified as Rho (D) positive or Rho (D) negative depending upon the presence or absence of D (Rho) antigen on them. D^u is the phenotypic term used to denote a weakened expression of the D (Rho). Detection of D^u (weak/partial D's) individuals is important as they can cause alloimmunization in Rho (D) negative individuals.

REAGENTS

Anti-A (Clone 11H5) and Anti-B (Clone 6F9) are ready to use solutions of the respective specific antibodies of the immunoglobulin class IgM prepared from the corresponding supernatants of mouse hybridoma cell cultures. Anti-D (Rho) (IgM + IgG) (Clone P3 x 61+MCAD6) is a ready to use solution containing specific antibodies of the immunoglobulin class IgM and IgG prepared from the supernatants of mouse hybridoma

cell cultures. These antibodies are a mixture of several monoclonal antibodies of the same specificity but having the capability of recognising different epitopes of the human red blood cell antigen A, B and D (Rho) respectively. Anti-D (Rho) (IgM+ IgG) does not detect the D^u's in the routine slide and tube test. For the detection of D^u's all samples negative with Anti D (Rho) (IgM+IgG) in the slide and tube tests should be further investigated by performing the D^u test procedure using the Anti-Human Globulin reagent.

Each batch of reagent undergoes rigorous quality control at various stages of manufacture for its specificity, avidity and performance.

REAGENT STORAGE AND STABILITY

1. Store the reagent at 2-8°C. DO NOT FREEZE.
2. The shelf life of the reagent is as per the expiry date mentioned on the reagent vial label.

PRINCIPLE

Human red blood cells possessing the corresponding antigen will agglutinate in the presence of antibody directed towards the antigen. Agglutination of red blood cells with Anti-A, Anti-B and Anti-D (Rho) (IgM + IgG) reagent is a positive test result and indicates the presence of the corresponding antigen. Absence of agglutination of red blood cells with Anti-A and Anti-B reagent is a negative test result and indicates the absence of the corresponding antigen. For D (Rho) and D^u (weak/partial D's) antigen, absence of agglutination of red blood cells by the D^u test procedure using Anti-D (Rho) (IgM+ IgG) and Anti-human globulin (Coombs) reagent indicates the absence of the corresponding antigen.

NOTE

1. In vitro diagnostic reagent for laboratory and professional use only. Not for medicinal use.
2. The reagent contains sodium azide 0.1% as preservative. Avoid contact with skin and mucosa. On disposal flush with large quantities of water.
3. Extreme turbidity may indicate microbial contamination or denaturation of protein due to thermal damage. Such reagents should be discarded.
4. Reagents are not from human sources, hence contamination due to HBsAg, and HIV is practically excluded.

SAMPLE COLLECTION AND PREPARATION

No special preparation of the patient is required prior to sample collection by approved techniques. Samples should be stored at 2-8°C if not tested immediately. Do not use haemolysed samples.

Anticoagulated blood using various anticoagulants should be tested within the below mentioned time period:

EDTA or HEPARIN	: 2 days
Sodium citrate or sodium oxalate	: 14 days
ACD or CPD	: 28 days

ADDITIONAL MATERIAL REQUIRED FOR SLIDE AND TUBE TESTS

Glass slides (60 x 85 mm), Test tubes (12 x 75 mm), Pasteur pipettes, Isotonic saline, Centrifuge, Timer, Mixing sticks, SBio Anti-Human Globulin reagent.

TEST PROCEDURE

Bring reagents and samples to room temperature before testing.

Slide Test

1. Place one drop each of Anti-A, Anti-B and Anti-D (Rho) (IgM + IgG) reagent on a clean glass slide.
2. To each reagent drop, add one small drop (50 µl) of whole blood.
3. Mix well with a mixing stick uniformly over an area of approximately 2.5 cm².
4. Rock the slide gently, back and forth.
5. Observe for agglutination macroscopically at two minutes.

Tube test

1. Prepare a 2-3% suspension of the red cells to be tested in isotonic saline.
2. Place one drop each of Anti-A, Anti-B, Anti-D (Rho) (IgM+IgG) reagent into corresponding labeled test tubes.
3. Pipette into each of the test tubes, one drop (~50 µl) of the test red cell suspension and mix well.
4. Centrifuge for 1 minute at 1000 RPM (125 g) or 20 seconds at 3400 RPM (1000 g) or alternative speed and time validated by the laboratory.
5. Gently resuspend the cell button, observing for agglutination macroscopically.

Microplate Method

SBio ABD Combipack reagents are standardized for use in Microplate technique, however it is recommended that each laboratory should standardize and validate their own procedure using "U" bottom microplates. The below mentioned procedure is a guideline which should be considered while standardization and validation process of procedure because of the variation in methods and equipment used in various laboratories.

1. Prepare 2-3% suspension of the red cells to be tested in isotonic saline.
2. Add 1 volume of Anti-A, Anti-B and Anti-D (Rho) (IgM+IgG) reagents using the reagent vial dropper into appropriate microwells.
3. Add equal volume of test red cell suspension.

4. Mix the contents of the microplate, taking care to avoid cross-well contamination.
5. Centrifuge the microplate at low spin (400g for 30 secs). Or Centrifuge at a speed and time as per the standardization and validation process done using same centrifuge.
6. Re-suspend the red cells using a microplate shaker or manually.
7. Read and record results.

Note: Incubation of 10-15 mins at Room Temperature may be included to enhance the reactivity.

D^uTEST PROCEDURE

1. Proceed with further testing of the Anti-D (Rho) (IgM+IgG) tube after step No.5 of Tube test.
2. Incubate the Anti-D (Rho) (IgM+IgG) tube at 37°C for 15 minutes.
4. Wash the contents of the tube thoroughly, atleast three times, with isotonic saline and decant completely after the last wash.
5. Add two drops of SBio Anti-Human Globulin (Coomb's) reagent and mix well.
6. Centrifuge for 1 minute at 1000 RPM (125 g) or 20 seconds at 3400 RPM (1000 g).
7. Very gently, resuspend the cell button, observing for agglutination macroscopically.

INTERPRETATION OF RESULTS

Slide and tube tests

Agglutination is a positive test result and indicates the presence of the respective red blood cell antigens. Do not interpret peripheral drying or fibrin strands as agglutination. No agglutination is a negative test result and indicates the absence of red blood cell antigens.

D^uTest Procedure

- (a) Agglutination in the D^u test indicates presence of the D^u antigen (weak / partial D's). No agglutination in the D^u test indicates absence of the D^u antigen (weak / partial D's). Negative reactions obtained in D^u test must be validated:- add 50µl of Coomb's control cells to the reaction mixture. A positive reaction confirms the activity of the Coomb's reagent and validates the negative reaction before the addition of the Coomb's control cells.
- (b) Mixed field agglutination in the D^u test on red cells from a recently delivered woman may indicate a mixture of maternal Rho (D) negative and fetal Rho (D) positive blood.
- (c) Red cells demonstrating a positive direct antiglobulin test cannot be accurately tested for D^u antigen (weak / partial D's).

NOTE

1. It is strongly recommended that red cells with known ABO & Rh characteristics should be periodically run, preferably on a daily basis to validate the reagent performance.
2. The quality control of reagents is performed using immediate spin method for Anti-A, Anti-B, Anti-A,B reagents.

3. The quality control reagents is performed by incubating at 37°C to 30 minute incubation for Anti-D IgM, Anti-D IgM +IgG reagents.
4. The reagents are expected to exceed minimum specifications/ acceptance criteria as per label claim for titre, specificity & avidity as laid down by transfusion medicine technical manual-2003 with reference to NIBSC-UK & WHO international; standards.

REMARKS

1. (a) Anti-A reagent react with Tn. Tn positive person must be excluded from donating blood as the occurrence of Tn is considered to be a symptom of a preleukaemic state and the red blood cells are polyagglutinable.
(b) Anti-B is truly negative reacting with acquired B characteristics.
2. In the tube test procedure, it is recommended that tubes with negative reactions should be recentrifuged and results read again after 5 minutes so that weak antigens are not overlooked.
3. As undercentrifugation or overcentrifugation could lead to erroneous results, it is recommended that each laboratory calibrate its own equipment and determine the time required for achieving the desired results.
4. Results of forward grouping obtained by using Anti-A and Anti-B reagents should always be reconfirmed by performing reverse grouping with known red cells.
5. It is strongly recommended that red cells with known ABO characteristics and known Rho (D) characteristics such as Rho (D) positive and Rho (D) negative red cells should be occasionally run, preferably on a daily basis so as to control reagent performance and validate test results.
6. After usage the reagents should be immediately recapped and replaced to 2-8°C storage.

7. The label minimum titre claim is based on A₁ group cells for Anti-A reagent, B group cells for Anti-B reagent and Rho (D) positive red cells (phenotype : R,r) for Anti-D (Rho)(IgM+IgG) reagent. This is based on the titration procedure as recommended by the manufacturer. Any deviation in test procedure could result in variable results.

PERFORMANCE CHARACTERISTICS

The performance of Anti-A, Anti-B and Anti-D (Rho) (IgM+IgG) complies with the common technical specifications of In-vitro diagnostic medical devices under the recommended methods. The performance of Anti-A, Anti-B and Anti-D (Rho) (IgM+IgG) was evaluated on over 3275 samples (from donors, clinical and neonates) drawn in the recommended anticoagulants. The evaluation demonstrated 100% specificity of each reagent and 100% sensitivity of Anti-A and Anti-B reagent. The sensitivity of Anti-D (Rho) (IgM+IgG) reagent by slide test is 99% and by tube test it is 99.53% of reagent versus the expected results with common known Rhesus phenotypes.

WARRANTY

This product is designed to perform as described on the label and the package insert. The manufacturer disclaims any implied warranty of use and sale for any other purpose.

BIBLIOGRAPHY

(1) Kohler C. & Milstein C. (1975), Continuous cultures of fused cells secreting antibody of predefined specificity. *Nature*, 256, 495-497. (2) Lee H.H., Rouger P., Germain C., Muller A & Salmon C. (1983). The production and standardisation of monoclonal antibodies as AB blood group typing reagents. Symposium of International Association of Biological Standardisation on monoclonal antibodies. (3) Human Blood Groups, by Geoff Daniels, 1st Ed., Blackwell Science, Oxford 1995. (4) HMSO, Guidelines for Blood Transfusion Services, 2nd Ed., 1994.



Manufactured by:

TULIP DIAGNOSTICS (P) LTD.

Plot Nos. 92/96, Phase II C, Verna Industrial Estate,
Verna, Goa - 403 722, INDIA.

Regd. Office: Gitanjali, Tulip Block, Dr. Antonio Do Rego Bagh,
Alto Santacruz, Bambolim Complex P.O., Goa - 403 202, INDIA.
Website: www.tulipgroup.com Email: sales@tulipgroup.com

Manufactured for:



Singapore
Biosciences PTE Ltd.

11 Yishun Street 51, #04-23, The Criterion,
Singapore 767971